



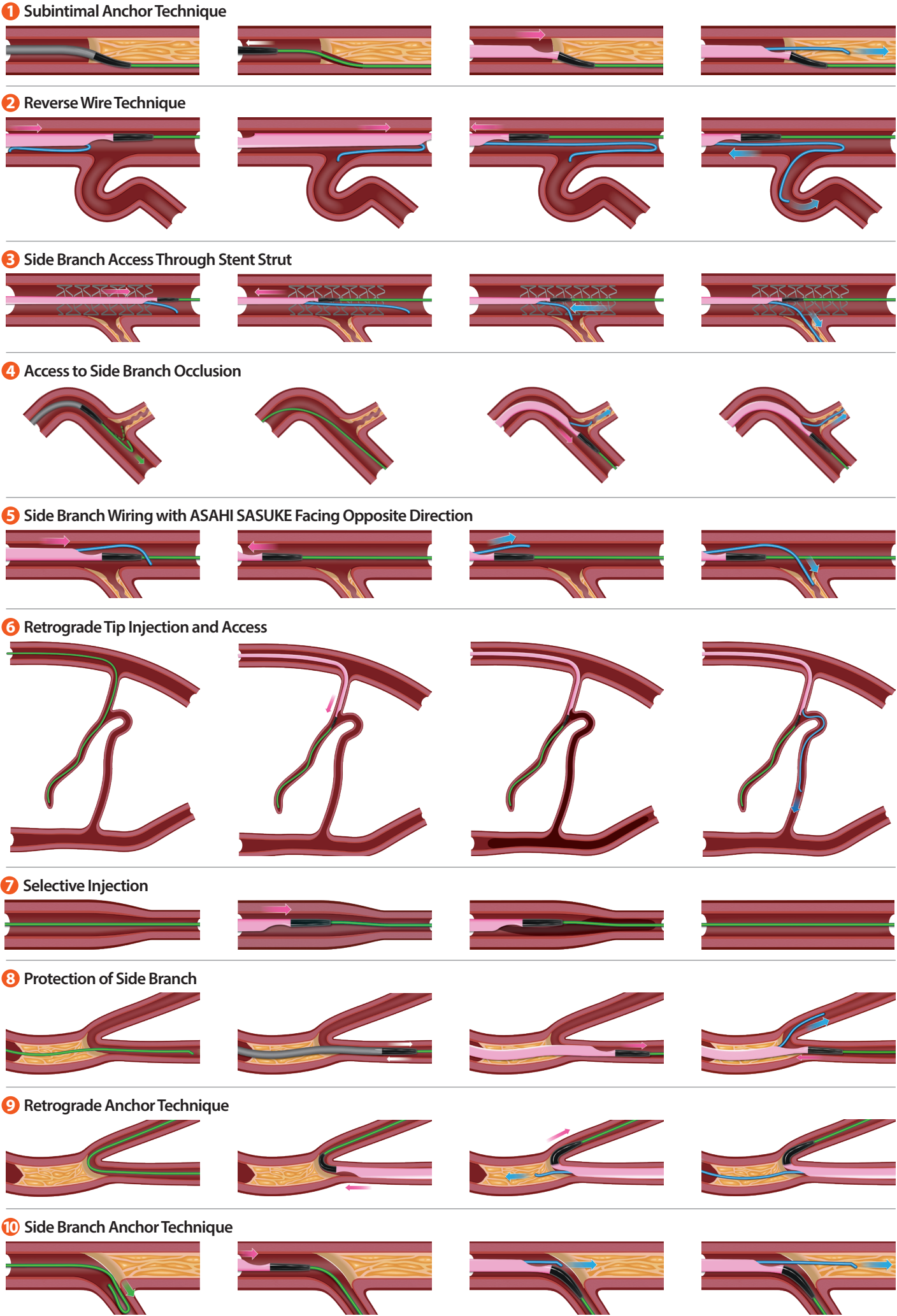

ASAHI SASUKE: the next generation of double lumen catheter

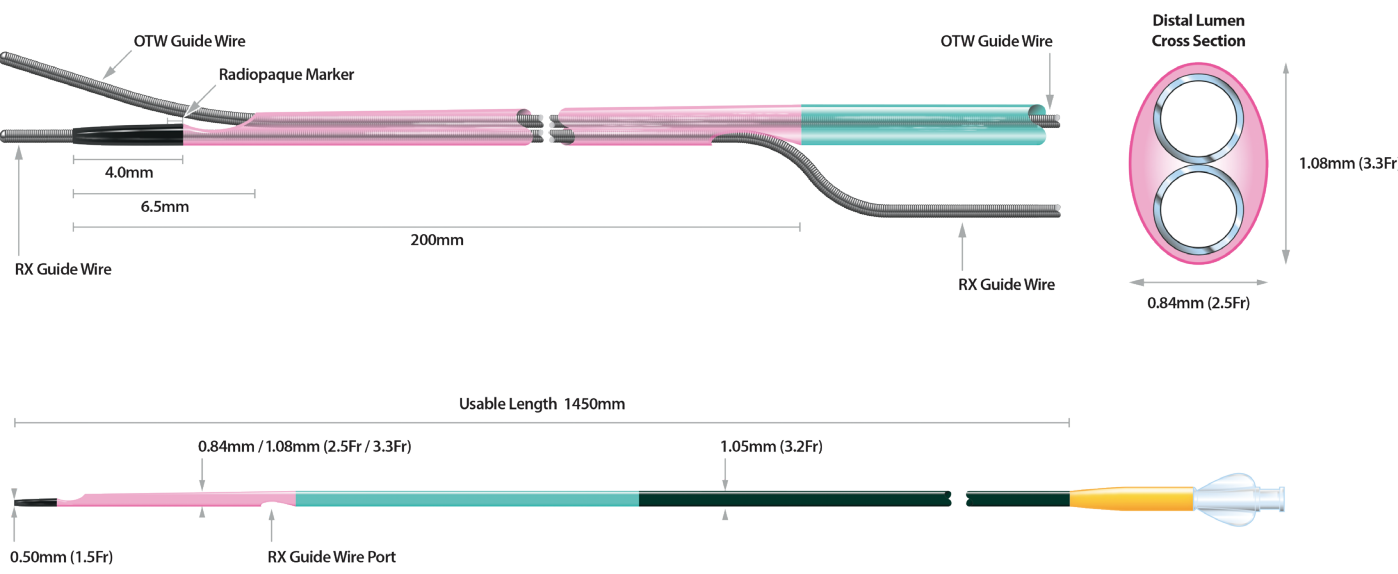

Ordering Information

Product	Catalogue No.	Outer Diameter			Inner Diameter		Usable Length	Recommended Guide Wire	Hydrophilic Coating Length
		Tip	Distal Shaft	Proximal Shaft	Tip	Shaft			
ASAHI SASUKE	SA145-33N	0.50mm (1.5F)	0.84mm/1.08mm (2.5F/3.3F)	1.05mm (3.2F)	0.40mm (0.016")	0.43mm (0.017")	145cm	0.36mm (0.014")	38cm



© 2018 ASAHI INTECC CO., LTD.
"ASAHI" and "SASUKE" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:



ASAHI INTECC CO.,LTD.
3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
Contact phone number: +81-561-86-9101
<http://www.asahi-intecc.com>

Distributed by:


ASAHI SASUKE/ver.2/L1/AMC-P18045






ASAHI SASUKE



Advantages of ASAHI SASUKE

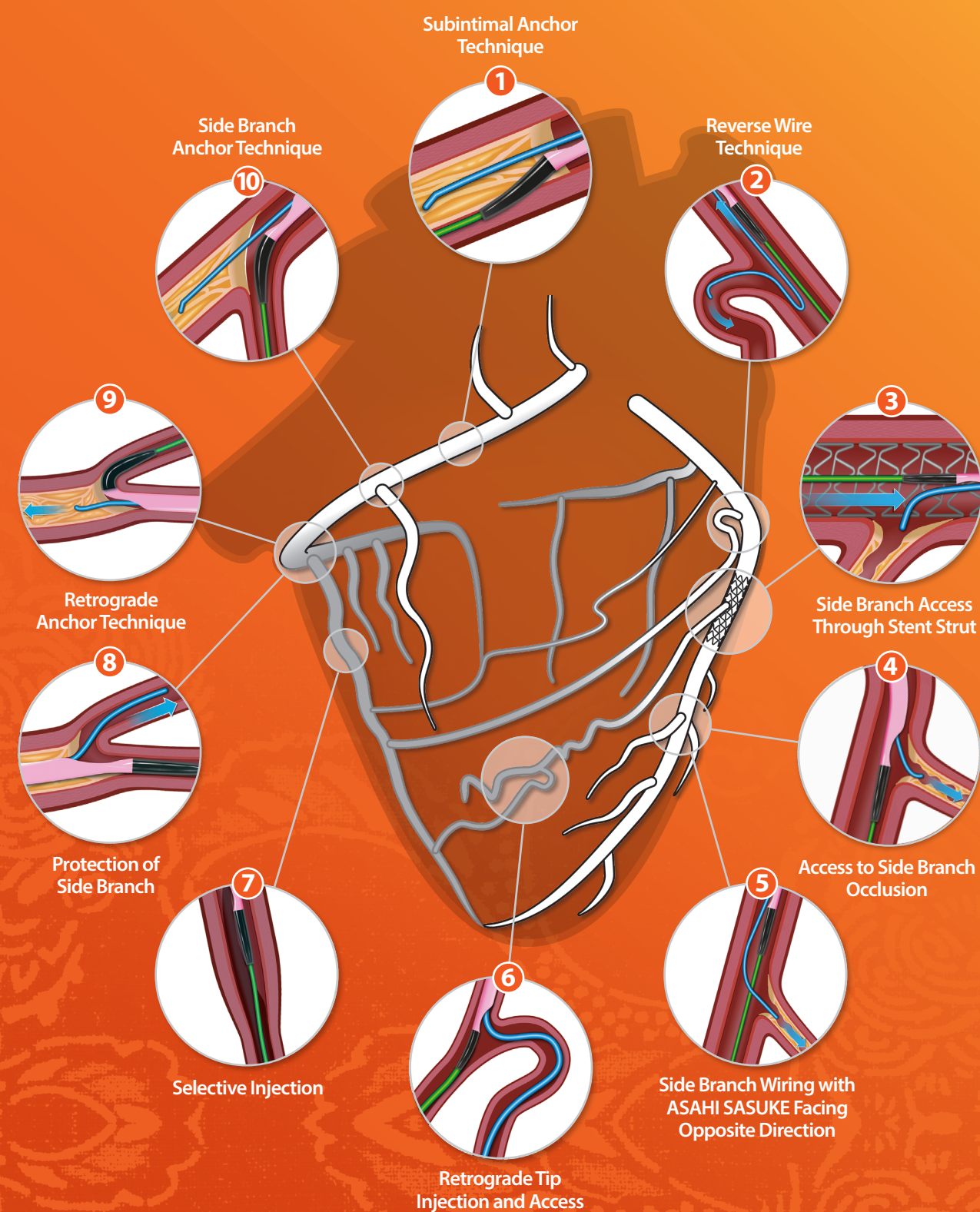

ASAHI SASUKE: the next generation of double lumen catheter

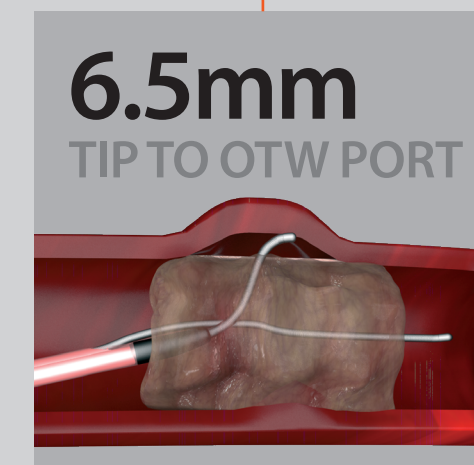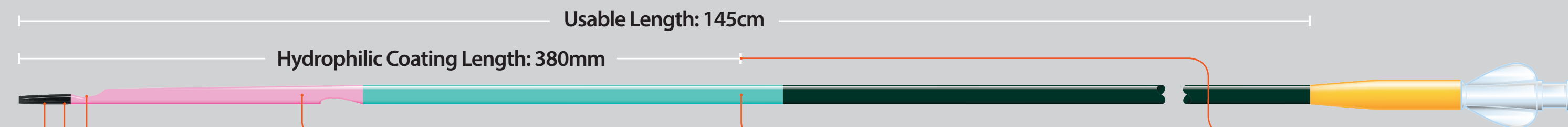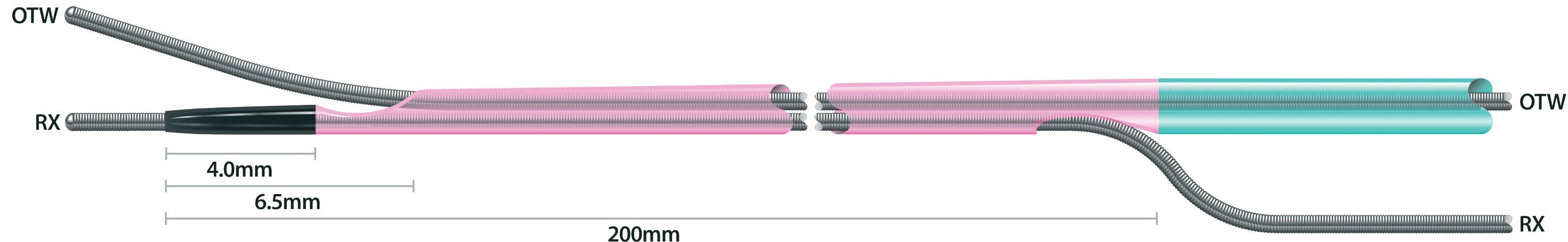


6.5mm tip to OTW port – negates the need for deep engagement inside branches and CTOs




Tapered soft tip - provides high trackability




Visible exit - know exactly where the wire is exiting

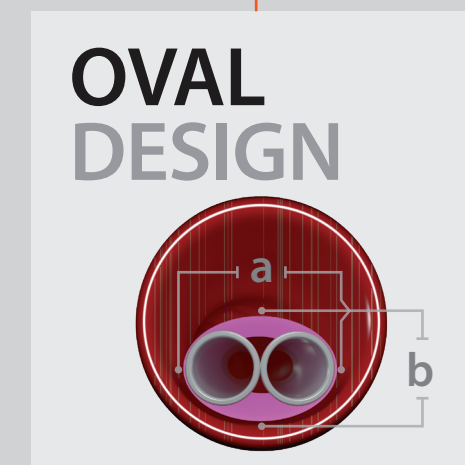

Oval design –
a) 1.08mm (3.3Fr) x
b) 0.84mm (2.5Fr)

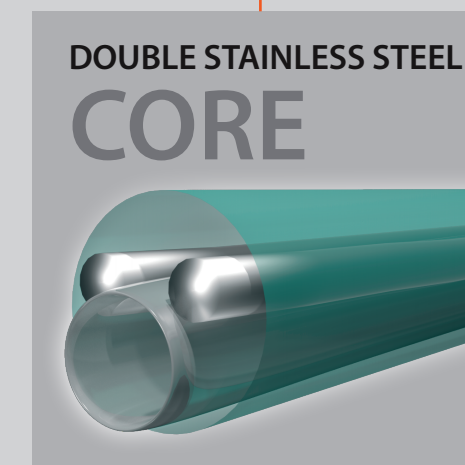

Double stainless steel core - kink resistance and pushability

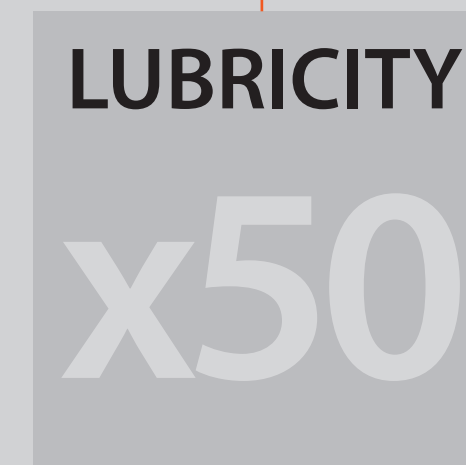

Lubricity - L³ Coating (long lasting lubricity) maintains trackability